



MEDIT *iE*

Start with Medit. **Evolve** with Medit.

iE, E for

Evolve

Easy Start

Scan your first case on installation day.

- Simple Plug & Scan
- Easy, intuitive UI/UX anyone can scan with
- No onboarding for existing Medit users

Excellence in Every Scan

Proven over 100M+ cases, consistent in every case.

- The same scan performance as Medit's flagship
- Consistent results in every case
- Performance that keeps improving with free updates

Expandable Workflow

Start with Scan & Send, expand within Medit when ready.

- Realtime Intelligence for scan and prep checks, when you need it
- Medit Apps for consultation and design as you grow
- Move up to the i900 Series · same account, same data

Ecosystem You Can Trust


With Medit Link, the cloud platform trusted by 200K+ users.

- HIPAA and GDPR regional data security and compliance
- True Open Platform · 150+ partner integrations
- A global service network you can rely on

Start with Medit. **Evolve** with Medit.



MEDIT *iE*

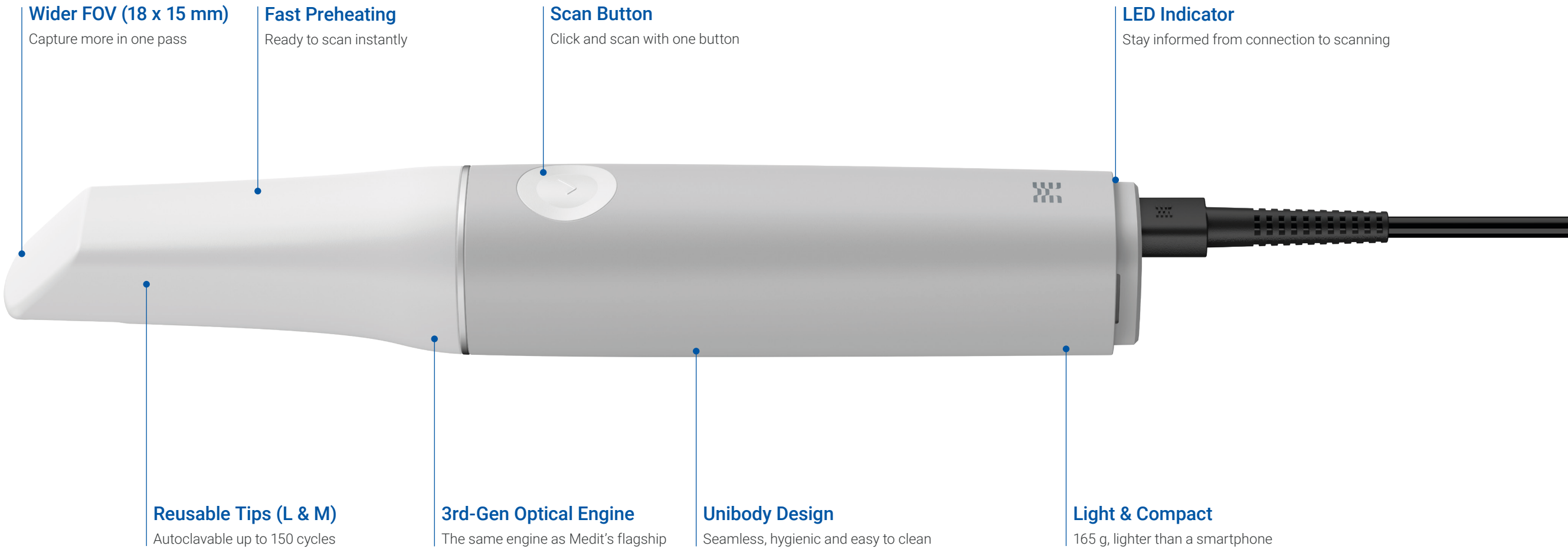


Optics & Scan

Sharper, wider, smarter. The optics behind every Medit scan.

3rd-Gen Optics, Sharper Than Ever

Advanced 10-bit imaging for ultra-realistic scan results.

- Enhanced clarity with higher bit depth
- True-to-life color reproduction
- High-resolution texture recognition

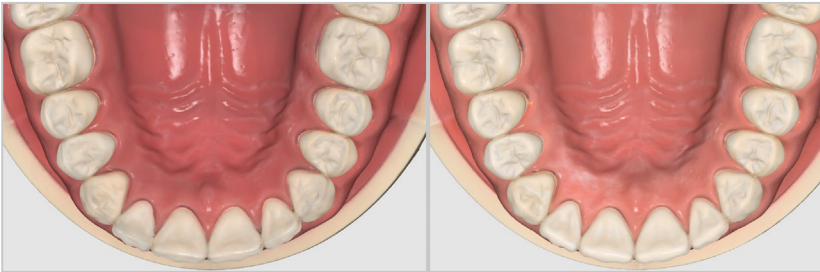


Medit 10-bit imaging technology

Comparative Product

Wider View, Faster Workflows

18 × 15 mm FOV and 30 mm depth for broader, deeper capture.

- Capture more per pass
- Ideal for complex or posterior regions
- Reduces scan time

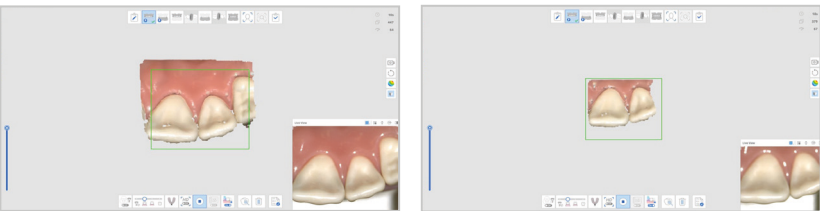


Medit iE

Comparative Product

Sharper, Smarter Metal Scanning


Built-in smart filters for accurate scans on reflective surfaces.

- Metal optimized filtering
- Reduced noise, clearer surfaces
- Reliable even in complex prosthetic cases

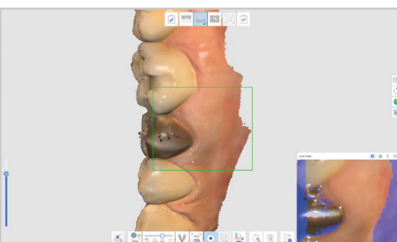

Scan & Send


From case creation to scanning to lab transfer.
Everything everyday dentistry needs, included from day one.

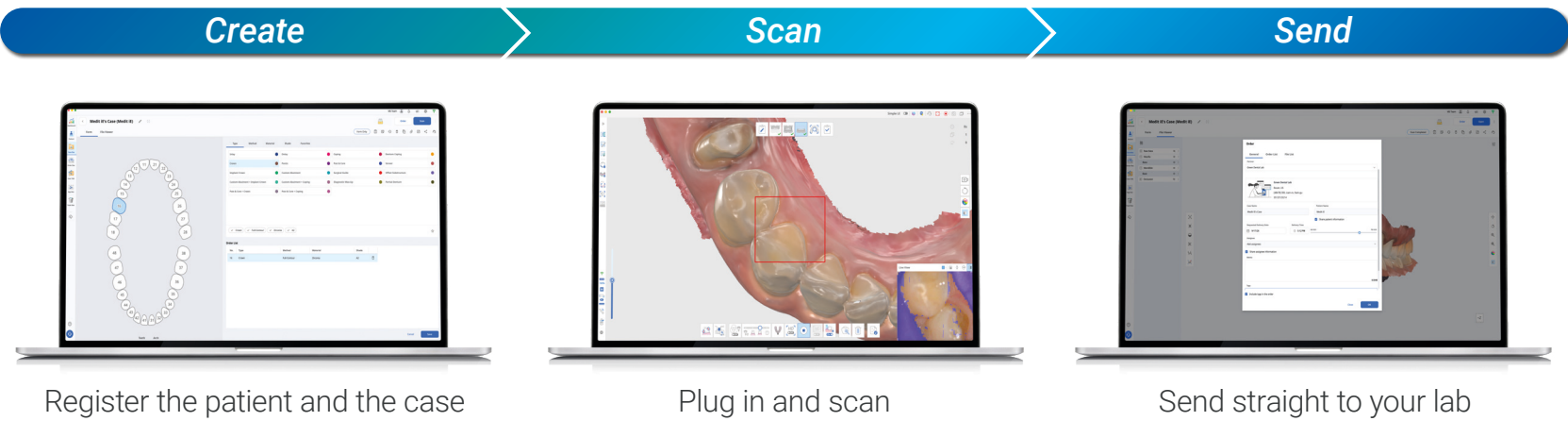

Medit Link

The cloud platform trusted by dental professionals worldwide

- Case management to lab communication, all in one place
- Data protected in line with HIPAA, GDPR and regional security standards
- Free software updates bring new features and better performance
- Works with 150+ partner software through open formats

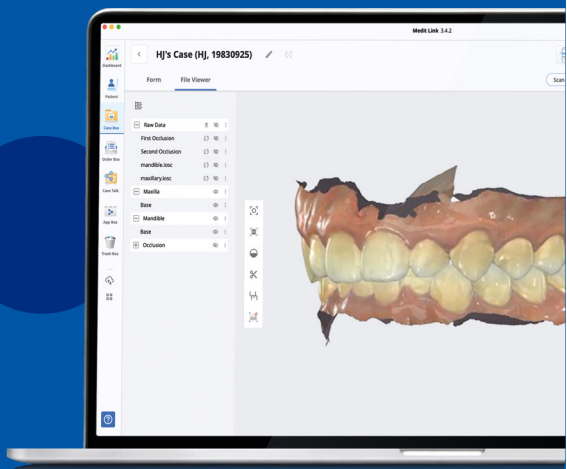

Expand Within Medit

Everything you add — and everything you move up to — stays in one Medit ecosystem.

Medit’s Latest Scan Technology, With the Right Configuration for Everyday Dentistry.

Add Software — Anytime, on This Scanner



Realtime Intelligence

Real-time scan & prep quality checks — while the patient is still in the chair.



Medit Apps

Patient consultation and in-house design — add only what your practice needs.

Added in Minutes. No New Hardware

**Optional add-ons, sold separately*



When You’re Ready for a Broader Range of Care

i900 Series

All-on-X Scanning

Full-arch implant workflows with Medit SmartX.

Medit Link Express

Scan wirelessly — patients follow along on iPad. Medit’s newest experience, effortless to start.



For Workflows Like These,
the i900 Series Is the Right Fit.

One Medit Account — Your Cases, Data and Lab Connections Move With You.

Medit iE Specifications

Category		Medit iE
Type		Wired
Imaging Technology		3D-in-motion video technology / 3D full color streaming capture
HD Scan		Yes
Light Source		LED
Scanning Frame		70 FPS (Maximum)
Anti-fogging Technology		Adaptive anti-fogging
Powder Required		No
Handpiece	Dimension	223.4 × 36.7 × 35.3 mm (L × W × H)
	Weight	165 g
	Dimension - Full Tip	Large: 36.1 × 34.1 × 90.8 mm (W × H × L) Medium: 36.1 × 34.1 × 90.4 mm (W × H × L)
	Dimension - Tip Head	Large: 26.9 × 19.7 mm (W × H) Medium: 22.4 × 16.3 mm (W × H)
Reusable Tip	Autoclavable	Up to 150 cycles
		Sterilize in a gravity autoclave at 135°C (275°F) for 10 minutes and dry for 30 minutes.
		Sterilize in a pre-vacuum autoclave at 134°C (273.2°F) for 4 minutes and dry for 20 minutes.
		Sterilize in a gravity autoclave at 121°C (249.8°F) for 30 minutes and dry for 15 minutes.
Connectivity	Interface	USB 3.1 Gen 1 (Power Delivery)
	Cable	2.0 m, detachable
Software		Medit Link (Windows / Mac)

Disclaimer

Product availability, regulatory approval status, and release timelines may vary by country. For accurate and up-to-date information, please contact your local authorized sales representative.

Why Medit

#1

Brand in user ratings*

Intuitive software, frequent free updates and an open ecosystem. The reasons users gave.

100M+

Cases scanned with Medit

Scan performance proven across more than 100 million cases.

200K+

Dental professionals

Clinics and labs on Medit Link, the cloud platform.

150+

Partner software integrations

A true open platform. Your data is never locked in.



190+

Countries



\$0

Software updates, always free



HIPAA · GDPR

Regional data security and compliance

The Brand That Opened Digital Dentistry to Everyone.

*InsightsZ, Global IOS Category Series, Sentiment Analysis of Leading IOS Brands Vol. I (2025)

